


Amnioram

Ago per amniocentesi
Amniocentesis needles



Per il prelievo transaddominale di liquido amniotico (eseguito in genere tra la 15ª e la 18ª settimana di gravidanza) per la diagnosi prenatale di difetti genetici del feto. Utilizzato anche per esami citogenetici e biochimici. Il mandrino è provvisto di punta eco-marcata che permette un inserimento ecoguidato. La cannula è provvista di un fermo scorrevole per determinare la profondità di inserimento dell'ago e di un raccordo luer lock per la raccolta del liquido amniotico a mezzo siringa. **AMNIORAM** è disponibile in 12 diverse misure.

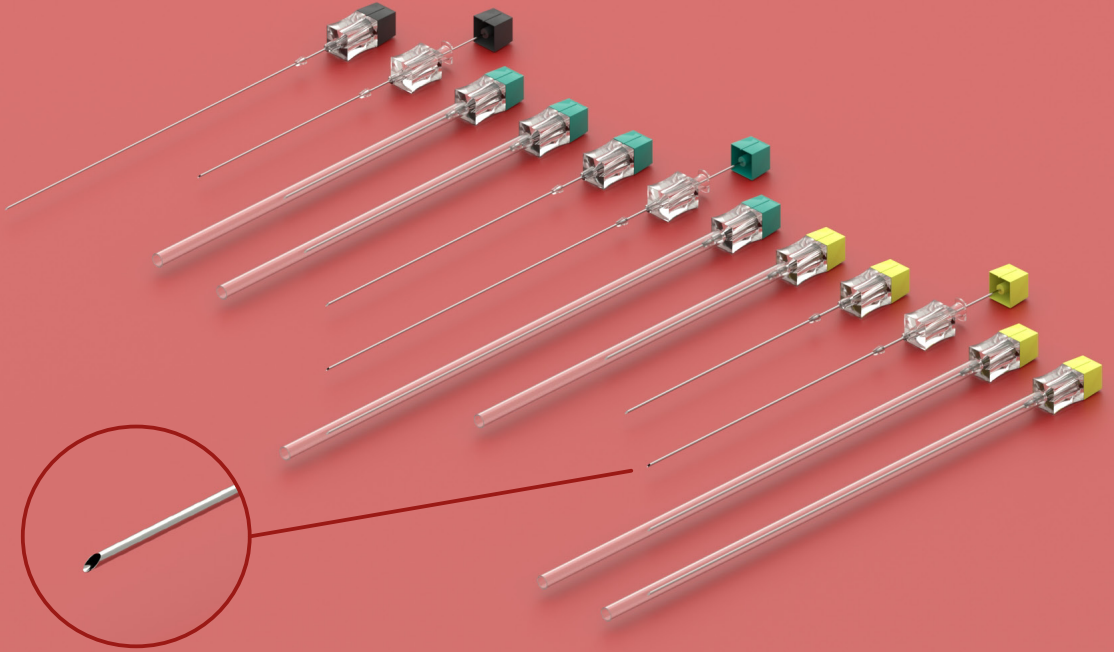

For transabdominal sampling of amniotic fluid (usually performed between the 15th and 18th pregnancy week) for prenatal diagnosis of genetic disorders of the foetus. It is also used for cytogenetic and biochemical studies. The mandrel is provided with an ultrasound detectable tip allowing an echo-guided insertion. The needle is provided with an adjustable stopper to determine insertion depth, and with a luer lock connector to collect the amniotic fluid through a syringe. **AMNIORAM** is available in 12 different sizes.

Cod.	Colore - Color	Descrizione / Description	Pz/Scatola - Pcs/Box
770212	Giallo / Yellow	20 G x 80 mm	25 pz./pcs.
770228	Giallo / Yellow	20 G x 100 mm	
770222	Giallo / Yellow	20 G x 120 mm	
770224	Giallo / Yellow	20 G x 150 mm	
770201	Giallo / Yellow	20 G x 180 mm	
770226	Verde / Green	21 G x 80 mm	
770227	Verde / Green	21 G x 100 mm	
770223	Verde / Green	21 G x 120 mm	
770202	Verde / Green	21 G x 150 mm	
770229	Verde / Green	21 G x 180 mm	
770233	Nero / Black	22 G x 90 mm	
770232	Nero / Black	22 G x 150 mm	

Fornito in confezione singola - monouso - sterilizzazione a EtO - marcato CE 0051

Supplied in single packaging - disposable - EtO sterilization - CE 0051 marked

