



Amniocent


Ago per amniocentesi con cannula centimetrata
Amniocentesis needles with centimeter-markings

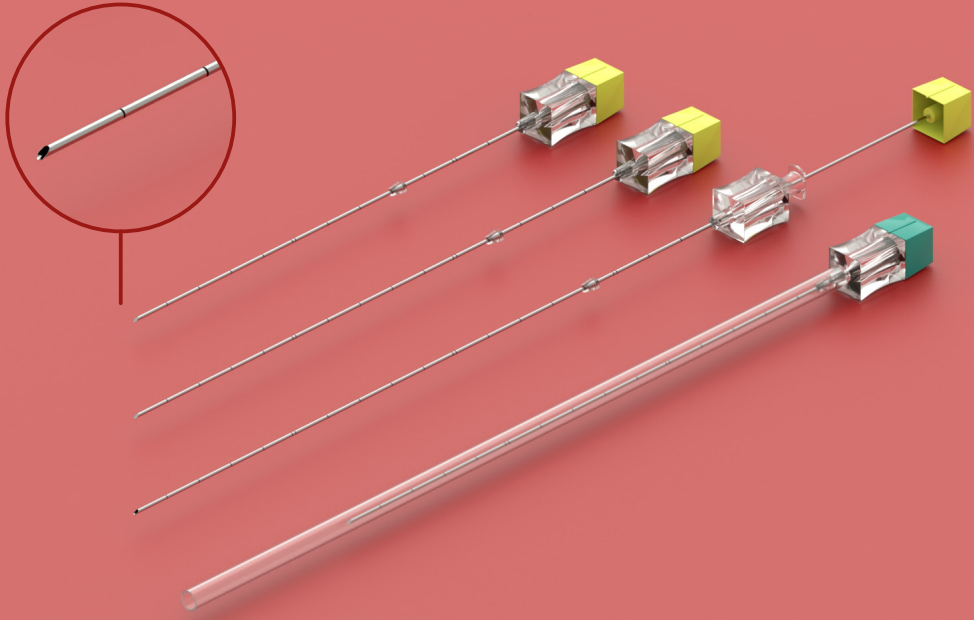

 Per il prelievo transaddominale di liquido amniotico (eseguito in genere tra la 15ª e la 18ª settimana di gravidanza) per la diagnosi prenatale di difetti genetici del feto. Utilizzato anche per esami citogenetici e biochimici. Il mandrino è provvisto di punta eco-marcata che permette un inserimento ecoguidato. La cannula è provvista di un fermo scorrevole per determinare la profondità di inserimento dell'ago e di un raccordo luer lock per la raccolta del liquido amniotico a mezzo siringa.

AMNIOCENT è disponibile in 4 diverse misure. La cannula è centimetrata (1 cm di risoluzione)* e presenta una doppia tacca di riferimento a 5 cm.


* Fornisce una misura approssimativa.

 For transabdominal sampling of amniotic fluid (usually performed between the 15th and 18th pregnancy week) for prenatal diagnosis of genetic disorders of the foetus. It is also used for cytogenetic and biochemical studies. The mandrel is provided with an ultrasound detectable tip allowing an echo-guided insertion. The needle is provided with an adjustable stopper to determine insertion depth, and with a luer lock connector to collect the amniotic fluid through a syringe.

AMNIOCENT is available in 4 different sizes. The needles is centimeter-marked (1 cm resolution)* and has a double notch at 5 cm.

* It gives an approximate measurement.

Cod	Colore - Color	Descrizione / Description	Pz/Scatola - Pcs/Box
772007	Giallo Yellow	20 G x 120 mm	25 pz./pcs.
772004	Giallo Yellow	20 G x 150 mm	
772006	Giallo Yellow	20 G x 180 mm	
772001	Verde Green	21 G x 150 mm	

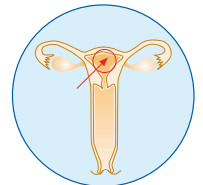

Fornito in confezione singola - monouso - sterilizzazione a EtO - marcato CE 0051

Supplied in single packaging - disposable - EtO sterilization - CE 0051 marked