



# IUD Extractor

Dispositivo per la rimozione della spirale  
Device for the removal of retained IUD

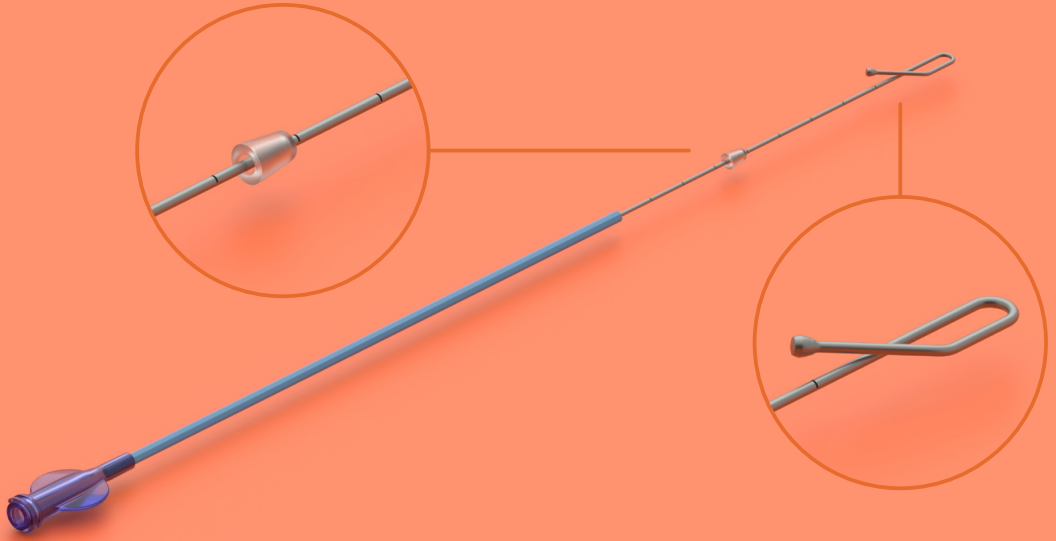


 Indicato per la rimozione di IUD nei casi in cui non sia visibile il filo di repere o per le spirali che ne sono prive (tipo anello cinese).

Il design della punta ne permette un facile inserimento e un agevole aggancio del filo o della struttura della spirale da rimuovere.

L'olivetta al termine della punta evita traumi tissutali.

La cannula in acciaio presenta una scala di misura di 9 cm (risoluzione 1 cm)\* ed è provvista di un fermo scorrevole per determinare la profondità di inserimento.

\* Fornisce una misura approssimativa.

 Designed for removing IUDs when the string is not visible or in presence of IUDs without strings (such as chinese ring).

The design of the tip makes the instrument easy to insert and allows simplified hooking of the IUD string or structure.

The ball at the end of the tip prevents tissue trauma.

The steel cannula has a 9 cm graduation marking (resolution 1 cm)\* and is provided with a sliding positioner to determine the insertion depth.

\* It gives an approximate measurement.

Cod.	Descrizione / Description	Pz/Scatola - Pcs/Box
770123	Lunghezza 28,5 cm Length 28.5 cm	10 pz./pcs.

Fornito in confezione singola - monouso - sterilizzazione a EtO - marcato CE 0051

Supplied in single packaging - disposable - EtO sterilization - CE 0051 marked

