







50 D // 80 D





SURTRON® 50 D is a radiofrequency electrosurgical unit suitable for minor monopolar surgery.

SURTRON® 50 D through the functions selection makes pure cut, cut coagulation (CUT 1 - CUT 2) and pure coagulation possible. The on-board microcontroller and digital display guarantees consistent safety and reliability under any conditions. SURTRON® 50 D offers outstanding professional surgical results due to its ergonomic and superior safety features.

Reference code Maximum output power CUT Maximum output power CUT/COAG1 Maximum output power CUT/COAG2 Maximum output power COAG Working frequency Patient circuit Selectable input voltage Mains frequency Electrical input power Size WxHxD mm Weight



SURTRON® 80 D is a radiofrequency electrosurgical unit suitable for minor monopolar surgery.

SURTRON® 80 D has a multi-function selector which makes cutting (CUT), cut coagulation (BLEND), coagulation (COAG) and micro-coagulation (MICRO) possible. The one-board microcontroller and digital display guarantees consistent safety and reliability under any conditions. SURTRON® 80 D offers outstanding professional surgical results due to its ergonomic and superior safety features.

Reference code 50 W – 400 Ω Maximum output power CUT Maximum output power BLEND 45 W - 400 Ω40 W – 400 Ω Maximum output power COAG 40 W – 400 Ω Maximum output power MICRO Working frequency Patient circuit 115/230 Vac Selectable input voltage Mains frequency Electrical input power 190x85x239 Size WxHxD mm Weight

10100.051

600 kHz

50-60 Hz

280 VA

2,5 Kgs

10100.081 80 W – 400 Ω  $60 \text{ W} - 400 \Omega$ 50 W - 400 Ω 20 W – 600 Ω 700 kHz 115-230 Vac 50-60 Hz 220 VA 190x85x239 2,5 Kgs







#### SAFETY

EN60601-1 EN60601-1-2 EN60601-2-2 Electrical Class: I CF MDD 93/42/EC Class: II b

#### CONTROLS

Patient/Plate circuit monitoring (80D) Output power monitoring Self check control

### 80 // 120 // 160



 $SURTRON^{\circ}$  80 / 120 / 160 is a high frequency electrosurgical equipment which is suitable to light and medium surgery.

SURTRON® 80 / 120 / 160, through its performance, allows pure CUT, cut-coagulation BLEND, superficial coagulation FORCED COAG, deep coagulation in absence of necrosis SOFT COAG and, with a special adapter, BIPOLAR coagulation.

The digital reading of delivered power and the monitoring of operative function by microcontroller guarantees the absolute reliability of working conditions.

 $\rm SURTRON^{8}~80~/~120~/~160~allows~a~highly~professional~surgery~thanks~to~the~user-friendly~and~safety~solutions~normally~used.$ 

The connection of neutral electrode is constantly monitored. Safety control of patient/plate contact using split neutral electrode. The possibility to control by the handle the monopolar output functions as well as the delivery of output power, allows to implement the surgical operation without turning away the surgeon attention from the surgical field.

Reference code
Maximum output power CUT
Maximum output power BLEND
Maximum output power FORCED COAG
Maximum output power SOFT COAG
Maximum output power BIPOLAR
Working frequency
Patient circuit
Selectable input voltage
Mains frequency
Electrical input power
Size WxHxD mm

Weight

10100.201
120 W-250 Ω
90 W-200 Ω
80 W-150 Ω
60 W-100 Ω
40 W-100 Ω
600 KHz
F
115-230 Vac
50-60 Hz
300 VA
254x104x288
5 Kgs

10100	.301
160 W	-250 Ω
120 W	-200 Ω
100 W	-150 Ω
80 W-	100 Ω
60 W-	100 Ω
600 KI	Ηz
F	
115-23	30 Vac
50-60	Hz
350 VA	A
254x1	04x288
5 Kas	



### FLASH // FLASH HF



SURTRON® FLASH HF is a High Frequency RadioSurgical equipment which is suitable to precision monopolar and bipolar surgery and micro-surgery without tissue alterations.

SURTRON® FLASH can supply RF operative powers either for programmable time as short as hundreds of a second.

SURTRON® FLASH through its performances allows pure CUT, cut-coagulation BLEND, incision with reduced prodution of eschar ENHANCED, superficial coagulation FORCED COAG, deep coagulation in absence of necrosis SOFT COAG, BIPOLAR coagulation and ABLATION. The digital reading of the delivered power and the overseeing through microcontroller of the operational functions, assure the absolute reliability of the conditions of job.

SURTRON® FLASH allows a highly professional surgery thanks to the user-friendly and safety solutions normally used.

The connection of neutral electrode is constantly monitored. Safety control of patient/plate contact using split neutral electrode. The possibility to control by the handle the output functions as well as the delivery of output power, allows to implement the surgical operation without turning away the surgeon attention from the surgical field.

## 120 160 HF 200

Reference code
Max output power CUT
Max output power BLEND
Max output power ENHANCED
Max output power FORCED COAG
Max output power SOFT COAG
Max output power BIPOLAR COAG
Max output power BIPOLAR ABLATION
Working frequency
Working frequency (1) FUNCTION ONLY
RF power pulse emission time
RF power pulse interval
Patient circuit
Selectable input voltage
Mains frequency
Electrical input power
Size WxHxD mm
Weight

10300.101
120 W-250 Ω
90 W-200 Ω
90 W-250 Ω
30 W-150 Ω
60 W-100 Ω
60 W-50 Ω
600 KHz
•
10 ms - 30 sec
1 - 9
=
115 - 230 Vac
50 - 60 Hz
300 VA
370x144x319
6 Kgs

10300.201
150 W-500 Ω <sup>(1)</sup>
100 W-500 Ω <sup>(1)</sup>
100 W-500 Ω <sup>(1)</sup>
100 W-150 Ω
80 W-100 Ω
60 W-50 Ω
60 W-100 Ω <sup>(1)</sup>
600 KHz
4.0 MHz
10 ms - 30 sec
1 - 9
F
F 115 - 230 Vac
115 - 230 Vac
115 - 230 Vac 50 - 60 Hz

10300.401
200 W-250 Ω
120 W-200 Ω
120 W-250 Ω
150 W-150 Ω
90 W-100 Ω
80 W-50 Ω
-
600 KHz
-
10 ms - 30 sec
1 - 9
F
115 - 230 Vac
50 - 60 Hz
350 VA
370x144x319
6 Kas

## 120 // // 200

#### **Standard Composition:**

- Reusable handle with switches 00205.00 - Handle for microsurgical needle 00201.02 - Assorted electrodes (10 pcs) 00500.00 Needles for micro-surgery (10 pcs) 500500.L11 Blade electrode 7 cm (3 pcs) 152-110 Needle electrode 7 cm (3 pcs) 152-120 Ball electrode 6 cm (3 pcs) Cable for neutral electrode Steel neutral electrode

Disposable split neutral electrode (2 pcs) Disposable neutral electrode (2 pcs)

Waterproof foot switch

Power supply cable 5mt

152-150 00404.08 5365 F7920 0350 00304.00

00100.01

AMBULATORY SURG







EN60601-1 EN60601-1-2 EN60601-2-2 Electrical Class: I CF MDD 93/42/EC Class: II b

#### CONTROLS

Patient/Plate circuit monitoring Output power monitoring Self check control



200



SURTRON® 200 is a high frequency electrosurgical equipment which is suitable to light and medium surgery.

 $SURTRON^{@}200\,allows\,to\,execute\,surgical\,minimally\,invasive\,procedure\,of\,resection, evaporation\,and\,coagulation.$ 

SURTRON® 200 through its performances allows pure CUT, cut-coagulation BLEND, incision with reduced prodution of eschar ENHANCED, superficial coagulation FORCED COAG, deep coagulation in absence of necrosis SOFT COAG and BIPOLAR coagulation.

The digital reading of the delivered power and the overseeing through microcontroller of the operational functions, assure the absolute reliability of the conditions of job.

SURTRON® 200 allows a highly professional surgery thanks to the user-friendly and safety solutions normally used. The connection of neutral electrode is constantly monitored. Safety control of patient/plate contact using split neutral electrode. The possibility to control by the handle the output functions as well as the delivery of output power, allows to implement the surgical operation without turning away the surgeon attention from the surgical field.

Reference code
Maximum output power CUT
Maximum output power ENHANCED
Maximum output power BLEND
Maximum output power FORCED COAG
Maximum output power SOFT COAG
Maximum output power BIPOLAR COAG
Working frequency

Working frequency
Patient circuit
Mains voltage
Mains frequency
Electrical input power
Size WxHxD mm
Weight

Patient to plate contact control
Split return electrode use allowed
Minimally invasive surgical treatment allowed
Monopolar and Bipolar indipendent output
Auto Start/Stop in Bipolar Coag

10100.401 200W - 250Ω 120W - 250Ω 120W - 200Ω 150W - 150Ω 90W - 100Ω 80W - 50Ω 600 kHz F 115-230 Vac 50-60 Hz 350 VA 370x144x319

6 Kgs

Working condition storing
Activation of the monopolar power with pedal and/or handle
Digital regulation and indication of the output power
Sound level control
Planning of the working condition

Minimally invasive surgery allowed



### **Standard Composition:**

Reusable handle with switches 00205.00 Assorted electrodes (10 pcs) 00500.00 Blade electrode 7 cm (3 pcs) 152-110 Needle electrode 7 cm (3 pcs) 152-120

Ball electrode 6 cm (3 pcs) Cable for neutral electrode 00404.08

Steel neutral electrode

Disposable split neutral electrode (2 pcs)

Waterproof foot switch

Power supply cable 5mt


152-150

5365

F7920 00304.00

00100.01



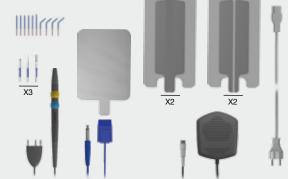

### CONTROLS

Patient/Plate circuit monitoring Output power monitoring Self check control

#### SAFETY

EN60601-1 EN60601-1-2 EN60601-2-2 Electrical Class: I CF MDD 93/42/EC Class: II b



### 00 HP // 400 F



#### TWIN OPERATIVE ACTIVITY

- Two output channels for monopolar cutting
- Two independent setting and output channels for the most frequent monopolar coagulation

#### **CONTROLLED CUTTING PROCEDURES**

- Automatic adaptive system of power\* (TACS - Tissue Adaptive Contact System)
- Enhanced cutting current for eschar free incision
- Cutting current with selectable coagulating property for bleeding controlled dissection
- · Low capacitive coupling by slow BLEND current for endoscopic procedures

#### TWIN CONTROLLED COAGULATIVE PROCEDURES

- SPEEDY\*\* coagulation and DEEP\*\*\*coagulation available to two operators
- SPRAY coagulation available by hand switch or foot switch controls

#### **WIDE BIPOLAR FEATURES**

- Bipolar cutting current with selectable bleeding control
- Slow BLEND cutting current
- Separately activable Automatic Start/Stop coagulation
- Vessel Sealing up to 7mm

#### **SURGEON AND PATIENT TOTAL SAFETY**

- Acceptability of contact impedance between patient body and return electrode continuously watched by Skin Plate Electronic Control (SPEC)
- H.F. Leakage currents continuously monitored
- Mains fluctuation automatically compensated before to start the delivery of HF power

#### PERSONALIZED PROCEDURES

- Independent choise of coagulative currents and power levels by each of two operators
- Ten personal configurable and immediately recallable procedures

#### **EXTENTION DRIVER AND INTERLOCK SIGNALS**

 Drived output connection to external Argon Plasma unit and Smoke Evacuator unit

#### **AUTOMATIC SELF DIAGNOSIS**

• Automatic diagnostic control of the actual complete functionality of the electronic unit and the connected accessories with coded report of failures

#### Reference code

Maximum output power monopolar pure CUT

Maximum output power monopolar cut-coag CUT80%

Maximum output power monopolar cut-coag CUT60%

Maximum output power monopolar cut-coag CUT40%

Maximum output power monopolar cut ENHANCED

Maximum output power monopolar cut-coag BLEND

Maximum output power monopolar coag SPEEDY Maximum output power monopolar coag DEEP

Maximum output power monopolar coag SPRAY

Maximum output power monopolar coag SPRAY ARGON

Maximum output power bipolar CUT

Maximum output power bipolar cut-coag CUT80%

Maximum output power bipolar cut-coag CUT60%

Maximum output power bipolar cut-coag CUT40% Maximum output power bipolar cut-coag BLEND

Maximum output power bipolar coag COAG

Maximum output power bipolar coag VESSEL SEALING

Working frequency Monopolar / Bipolar

Patient circuit

Selectable input voltage

Mains frequency

Electrical input power

Size WxHxD mm

Weight

470x150x400

17.5 Kas

0400.802	
00 W - 300 Ω	
50 W - 300 Ω	
00 W - 300 Ω	
50 W - 300 Ω	
00 W - 500 Ω	
00 W - 300 Ω	
00 W - 500 Ω	
00 W - 200 Ω	
0 W - 2000 Ω	
0 W - 2000 Ω	
20 W - 50 Ω	
20 W - 50 Ω	
20 W - 50 Ω	
00 W - 50 Ω	
20 W - 150 Ω	
20 W - 50 Ω	
00 W - 50 Ω	
25 kHz / 525 kHz	
15 - 230 Vac	
0 - 60 Hz	
000 VA	

400 W - 300 Ω 300 W - 300 Ω 250 W - 300 Ω 200 W - 300 Ω 250 W - 500 Ω 250 W - 300 Ω 120 W - 500 O 120 W - 200 Ω 100 W - 2000 Ω 100 W - 2000 Ω 120 W - 50 Ω 120 W - 50 Ω 120 W - 50 Ω 100 W - 50 Ω 120 W - 150 Ω 120 W - 50 Ω 200 W - 50 Ω 425 kHz / 525 kHz 115 - 230 Vac 50 - 60 Hz 1000 VA

470x150x400

17.5 Kas

<sup>\*</sup>Monopolar CUT mode until 100 Watt.

<sup>\*\*</sup>SPEEDY coagulation also said Fulgurate or Forced coagulation
\*\*\*DEEP coagulation also said Pin Point, Dessicate or Soft coagulation

AMBULATORY SURG SAFETY CONTROLS

EN60601-1 EN60601-1-2 EN60601-2-2 Electrical Class: I CF MDD 93/42/EC Class: II b Patient/Plate circuit monitoring Output power monitoring Self check control H.F. Leakage currents control Compensation of mains fluctuation

## 300 HP // 400 HP

#### **Standard Composition:**

- Disposable handle with switches (5pcs)
- Reusable handle with switches
- Blade electrode 7 cm (3 pcs)
- Blade electrode 16 cm (3 pcs)
- Needle electrode 7 cm (3 pcs)
- Ball electrode 6 cm (3 pcs)
- Cable for neutral electrode
- Disposable split neutral electrode (5pcs)
- Electrode cleaning sponge
- Double water-proof foot-switch (2 pcs)
- Power supply cable 5 mt


755VL 00205.00 152-110 152-115

152-120 152-150

00404.08 F7920

F7520 00305.04

00100.01

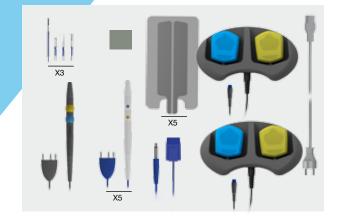



# SURTRON EVAC



**SURTRON EVAC** is an aspiration and filtering system for the smokes produced during the surgical procedures.

The smokes, generated during the surgical procedures, have a bad, strong and long smell. These smokes consist of water vapour, organic gas, visible or invisible solid particles and viral parts. So a good aspiration and filtration eliminates bad smells, reduce the bacteriologist and/or viral risks, and solves, the particulary annoying problem of the optimal vision of the zone of the intervention, during the laparoscopic procedures.

In **SURTRON EVAC** the activation of the aspiration can be:

AUTOMATIC, through electronic system of survey of the activation of the high frequency electrosurgical unit;

MANUAL, acting directly on the unit through a foot switch (optional)

The flux of the aspiration can be regulated to be suitable with the singular necessities. For laparoscopic procedures there is a specific function, quickly selectable.

**SURTRON** EVAC has an electronic control of the aspiration with the indication of the usury status of the filters.



Reference codes
Input voltage
Mains frequency
Electrical input power
Max vacuum flow rate
Max static suction
Filter type
Filter efficiency
Particles diameter
Size WxHxD mm
Weight

10200.10  $_{(230\ Vac)}$  - 10200.10/115V  $_{(115\ Vac)}$  230 Vac or 115 Vac 50 Hz 800 VA 1000 LPM (35 CFM) 250 mbar (83 in H<sub>2</sub>O) ULPA with active carbon 99,999X% 0,3 micron 370x144x319 4 Kg

- AUTOMATIC AND / OR MANUAL FUNCTIONING
- REGULATION OF THE ASPIRATION FLUX
- AUTOMATIC FUNCTION
  FOR LAPAROSCOPY
- CONTROL AND INDICATION OF USURY STATUS OF THE FILTER
- FAST SUBSTITUTION OF THE FILTERS
- DIFFERENT TYPE OF FILTERS
- ACTIVATED CARBONS FILTERS FOR THE ELIMINATION OF THE SMELLS



There are two different kind of filter, quickly substituible and constitued by a filtering ULPA type with a active carbons part dedicated to the elimination of the bad smells of the electrosurgical smokes and connector of 7/8" (22mm).

The filters have an electronic system for the control of the status of use.

### **FILTER BLUE**

The BLUE filter is the standard filter and it is included with the unit for a first time of use.

### FILTER RED

The **RED** filter is an optional spare part, created for a long time of use.

#### **Standard Composition**

- Power cable
- First ULPA Filter (5h) BLU
- Apiration kit: A Filter B Pipe dia. 22mm lenght 2,4 mt (100")
- 4 Handle apsiration kit: Adapter 10mm/22mm Pipe dia.10 Lenght 1,8 mt (72") Clip to fix the pipe on the handle

Adapter for the handle

#### **Replacement parts**

- ULPA Filter (20h) RED
- Aspiration kit (6 pcs)
- Handle aspiration kit (12 pcs)

00900.FU01 00900.02/06 00900.01/12



# **SURTRON® ABC**



**SURTRON** As a module that can be connected with a compatible electrosurgical unit, specific for coagulation electrosurgical procedures with Argon gas, indicated for both open surgery and laparoscopy.


Argon is an inert gas that used in combination with specific currents for high frequency coagulation generated by an electrosurgical unit, allows to obtain a superficial coagulation without contact thanks to the ionization of the gas flow through the high-frequency current.

**SURTRON** ABC is realized so that the connection to an electrosurgical unit with compatible characteristics to the Argon procedures, allows the erogation both through a manual activation and foot activation. The connection to the ABC unit doesn't preclude the use of other functions of the electrosurgical units, in fact, in case of use of the compatible LED units, the erogation of the Argon gas occurs only in case of erogation of current of high-frequency coagulation (SPRAY).



Reference code
Input voltage
Mains frequency
Electrical input power
Flow regulation
Maximum inlet pressure
Electrical Class
Size WxHxD mm
Weight

10500.20 90-240 Vac 50/60 Hz 70 VA 0.2 - 10 LPM 3.0 bar (43,5psi) I CF 470x150x400 8 Kgs

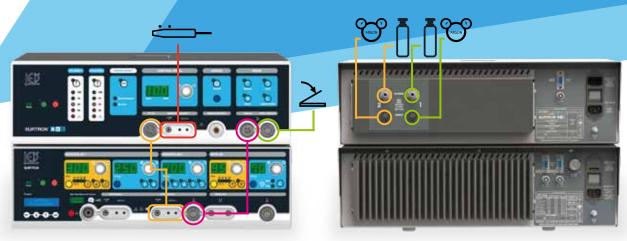


Connection with SURTRON HP (Not included)

- REGULATION AND CONTROL OF THE FLOW GAS SUPPLY
- STATUS CONTROL OF THE GAS ARGON CYLINDER
- ACTIVATION THROUGH THE HANDLE OR THE FOOT-SWITCH
- SPECIFIC FUNCTION FOR LAPAROSCOPY

- AUTOMATIC FUNCTION FOR COMPATIBLE



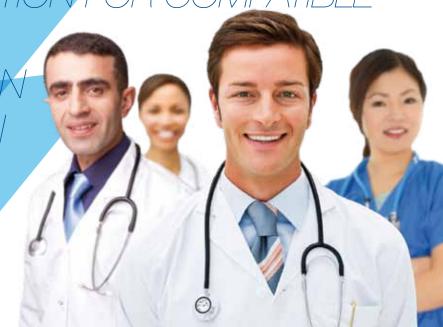

## www.surtron.com



LED SpA via Selciatella, 40 - 04011 - Aprilia (LT) - Italy www.led.it - info@led.it - www.surtron.com Tel. +39 06 92870045 Fax +39 06 92870046







